



shijiazhuang tiantian technology co.,ltd
hadermal-filler.com

HYALURONIC ACID
DERMAL FILLER
INJECTABLE HYALURONIC ACID

PDO Thread Lift

• A new thread lift technology with PDO thread lifting system
• A new thread lift technology with PDO thread lifting system

Microneedling Electric Dr. Pen Ultima M7 Micro Needling Derma Pen

Our Product Introduction

Basic Information

- Place of Origin: china
- Brand Name: dr pen
- Certification: ce
- Model Number: m7
- Minimum Order Quantity: 1
- Price: 55
- Packaging Details: one pen two needle
- Delivery Time: 3-5day
- Payment Terms: L/C, D/A, D/P, T/T, Western Union, MoneyGram
- Supply Ability: 10000




Product Specification

- Highlight: Electric Dr Pen Ultima M7 , Dr Pen Microneedling



More Images

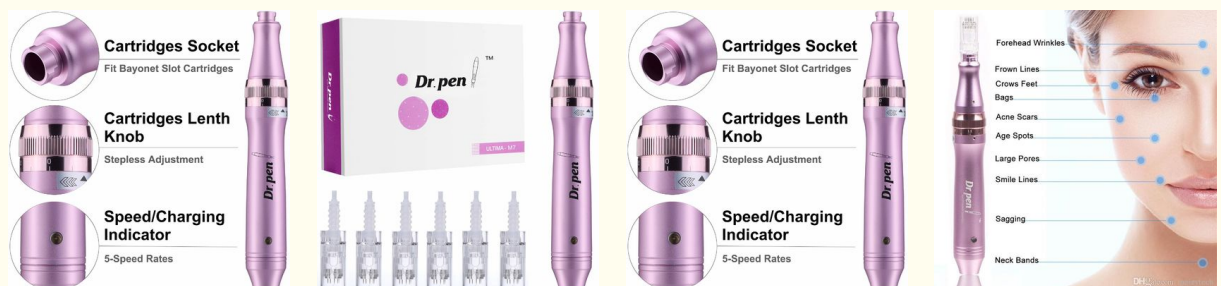

for more products please visit us on hadermal-filler.com

Product Description

Dr. Pen Ultima M7 Meso Micro Needling Machine Derma Pen

Dr Pen Ultima M7 microneedling pen is used mostly for face wrinkles, scars, stretch marks, and facial lines.

Dr Pen uses sterile, individually sealed, disposable needle cartridges for optimal safety and comfort.

The adjustable dial allows treatment that adapts to individual client's needs.

Dr Pen is more economical than the derma roller, you just need to recharge the pen and replace cartridges after every use.

Feature:

The AMTS (Auto Microneedle Therapy System) helps to:

Repair and revitalize complexion.

Reduce sagging skin, age spots, scarring (including acne scars), unattractive large pores & sun damage.

Reduce deep-set and fine-line wrinkles.

Restore skin's natural, youthful beauty.

Reduce cellulite & stretch marks.

Reverse age-related hair loss, restoring healthy, thick hair, a hallmark of youth.

Results rival those of aestheticians and cosmetic surgeons.

Mainly help to solve the following problems:

① Treatment of acne scars

② Scar repair

③ Smooth fine lines

④ Hair loss and hair treatment recovery

⑤ Stretch mark repair

⑥ Large pores

⑦ Tattoo removal

⑧ Reduction and removal of cellulite

⑨ Anti-edema

⑩ Hyperpigmentation & Melasma

Specification:

Product name: Dr pen Ultima M7

Working level: 5 levels of speed

Power type: Wired

Needle type: Bayonet needles

Application: For Beauty, Skin care, Home Use

plug-in model, must connect to the power when using

5 Gears Care Mode

Low Gear Activation

8000rpm/second

Open the skin channel and activate facial skin

Medium Gear Care

12000rpm/second

Layer-by-layer skin penetration, bring you infinite stimulation

High Gear Stimulation

16000rpm/second

Deep stimulation, essence of high efficiency infiltration



How to Operate:

1. Setting

Adapter Mode

(1) Plug in DC adaptor to the jack.

(2) Plug in DC adaptor to AC outlet.

2. Preparation

(1) Take out needle cartridge in the sterilized packing

(2) Fitting round slots 1 into 2 of the machine, push needle cartridge horizontally up to the end

(3) Lock the cartridge by turning it to the counter-clockwise

(4) Adjustment of needle length should be done after turning on the machine
(5) If you want to remove the needle cartridge, turn it to clockwise and push it out horizontally, fitting round slots (slot 1, slot 2) together
3. Adjustment of speed and needle length
(1) Adjustment of Speed and ON/OFF
(2) 5speed level control
Pressing the switch 2 seconds to open the device,
It starts from the lowest speed to highest speed when you turn on the power firstly
The speed cycles low to high. Please press the switch 2 seconds when you need to stop the device
The flicker frequency of high faster ,the speed is higher
(3) Adjustment of needles length
You can adjust the length of needle by turning the adjustment ring clockwise or counterclockwise
Turning it clockwise, length of it will be longer
Turning it counterclockwise, the length of it will be shorter

Attention:

For your Safety, use one roller per person only.

- 1.Do not force the roller into the skin but simply use mild pressure.
- 2.Roller obliquely at the same intensity covering the same areas vertically and horizontally.
- 3.Apply MTS protocols according to clients"s needs.
- 4.After the treatment, apply MTS post-care solutions for skin protection.
- 5.For safety and correct use, read the manual thoroughly before the roller.

Package includes:

- 1* Dr Pen Ultima M7 micro needling pen
- 2* 12-pin replacement cartridges
- 1* User Manual
- 1* Power adapter
- 1* storage box



shijiazhuang tiantian technology co.,ltd



008619931162089



info@tiantianbike.com



hadermal-filler.com

NO 260 HUANAN ROAD SHIJIAZHUANG CITY HEBEI CHINA