

## PRF Liquid Tube 10ml 15ml 13ml ACD Gel PRF Tube

Our Product Introduction

for more products please visit us on [hadermal-filler.com](http://hadermal-filler.com)

### Basic Information

- Place of Origin: original
- Brand Name: prp
- Certification: ce
- Model Number: 10ml 13ml 15ml
- Minimum Order Quantity: 1
- Price: 5
- Packaging Details: box
- Delivery Time: 2-3DAY
- Payment Terms: T, T/T, Western Union, MoneyGram
- Supply Ability: 10000



### Product Specification

- Name: PRP PRF Tube
- Use For: Anti Wrinkle
- Key Words: Tightening Gel
- Effecton: Butt Body Face Lift Remove Wrinkle
- Grade: Deep / Deep Plus
- Payment: TT/PayPal/Western Union/Credit Card
- Materials: 150 Mg Poly-L-lactic Acid (PLLA)
- Used For: Beauty Medical Clinic
- Keyword: Filorga 135
- Packing: Box
- Medical Type: Injection
- Dose: 150mg/5ml
- Pain Level: Mild Discomfort
- Recovery Time: 1-2 Days
- Volume: 5ml\*6 Vials



### More Images



### Product Description

PRP Tube 10ml with ACD Gel OEM for hair loss platelet Rich Plas Prp Centrifuge Tube With ACD Gel

Application	Skin rejuvenation, hair lose, dental, orthopedics
Vacuum available	8ml, 10ml, 12ml, 15ml, 20ml, 25ml, 30ml, 40ml, 50ml
PRP tube advantage	Sterile and pyrogen free
Brand Name	OEM PRP tube
Disinfecting Type	Irradiation
Material available	Glass or PET
Additive	Sodium citrate and Gel; ACD and gel
Shipping	DHL, TNT, UPS, EMS, Fedex, Air cargo and so on
Package	Packing in medical grade sterile blister

Our Product Introduction



#### SAFETY HEAD COVER

- Choosing high-quality synthetic butyl rubber stopper, which is very suitable for puncturing and sampling of automatic analyzer. It has the advantages of good air tightness, small puncturing force, few fallen scraps, not plugging holes and flicking needle.
- Patent design, butyl rubber stopper closely matches with safe cover so that it does not roll over and separate and the cap body is easy to open and close.
- Suitable for the use of various decapping centrifuges and cap opening machines.
- The colors of head cover conform to internationally accepted standards.

#### HIGH QUALITY TUBING

- Adopting high quality PET material with stable properties and good air tightness.
  - Inner wall of PET tube is treated with siliconization. It can avoid cell wall hanging and has extremely smooth surface.



#### HIGH EFFICIENT SEPARATION GEL

- Using high quality inert separation glue, has no physical and chemical interference with blood.
- Can get clear and transparent serum specimens after quickly centrifuging.
- High temperature resistance (45 °C); Convenient for freeze and storage; Has the properties of stability; rarely appear the phenomenon of "oil dropping" after centrifugation.
- Can be sterilized by  $\gamma$  irradiation.

#### HIGH QUALITY ADDITIVES


- The additives are all imported to avoid different batches.
- Complete variety and multiple forms; The spray drying technology makes anticoagulation uniform, soft and thorough.
- Presetting vacuum accurately to ensure that the proportion of blood sample and additive is accurate.



#### CUSTOMIZED LABELS

- Depending on customers' requirements, different materials and labels with specific logos can be customized. Printing prefabricated barcode is easy to identify and wear-resistant.
- Tearing type labels and barcode tag tubes have been stuck with customized labels before leaving the factory.

#### • Recommended Blood Collection Sequence:

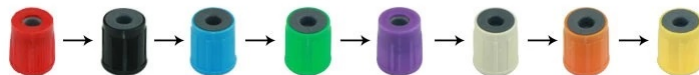

shijiazhuang tiantian technology co.,ltd



008619931162089



info@tiantianbike.com



hadermal-filler.com

NO 260 HUANAN ROAD SHIJIAZHANG CITY HEBEI CHINA